FR AB ENERGIZER MODE D'EMPLOI

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EN AB ENERGIZER INSTRUCTION MANUAL

Congratulations! You are now the proud owner of the most advanced electronic stimulation training unit currently available on the market. IMPORTANT: first read these instructions for use carefully before starting to use the Ab Energizer.

Ab Energizer Electronic Stimulation Units (ESU) issue soft electronic impulses in combination with the training belt that allow your muscles to contract and relax. Different muscles demand different degrees of contraction. It is, therefore, important to understand that the period of rest between the contractions (= the relaxation of the muscles) is the main part of the total workout. Our Ab Energizer is not a like other muscle stimulators. This advanced device uses two anatomically situated stimulators that will contract and relax in a mild and soft manner. You can really see your muscles move. Do not expect any vibrations or noise from the stimulation device when it is operational. The Ab Energizer is a mobile device that you can use anytime anywhere. The Ab Energizer is not a medical or therapeutic device.

How it works

Rub a generous amount of contact gel on the places that you want to train and attach the belt. You can use the Ab Energizer around your stomach, waist, thighs or buttocks. Attach the belt to the place where you rubbed in the gel and set the program and intensity of the ESU at the place that you will be training. The Ab Energizer will immediately focus on the relevant muscle group and will, thus, help to isolate and train the muscle group.

NOTE: DO NOT REMOVE LABELS FROM THE INSIDE OF THE TRAINING BELT.
IMPORTANT: BEFORE USING THE DEVICE PLEASE CONSULT YOUR DOCTOR.
NOT SUITABLE FOR CHILDREN.

You must press the units on the belt correctly and, next, position the belt with contact gel on your body. By starting the program, the Ab Energizer will complete its full program. We would, therefore, like to make you aware that the contact gel must always be generously rubbed on the places where you use the Ab Energizer for safety reasons. If little or no contact gel is used, it is possible that areas will appear that have been affected by friction. It is, therefore, advisable to apply the gel generously.

Mode Descriptions Indicator Light and Control Button

- Four Indicator Light: L1, L2, L3, L4
- ON/INC Button can turn off the power and gradually increase the intensity of vibration
- OFF/DEC Button can turn off the power and gradually reduce the intensity of vibration
- PROGRAM Button can switch between the six modes (Mode A/B/C/D/E/F)

Operating Instructions


1. With a CR2032 button battery.
2. Get a battery into the battery compartment (positive terminal down), close the battery door, open the plastic sheet on the viscous sticker, and then affix the massage patch to your desired position.
3. Press ON/INC button once and the indicator light on, then press PROGRAM button, so that select your desired working mode. Press the ON/INC again to increase intensity. Shut down or reduce the intensity, press the OFF/DEC button repeatedly.

Control Specifications


1. With 10 massage intensity and 6 massage working mode options.
2. LED display.
3. Normal working for about 15 minutes, then automatic power off
4. Press ON / INC button, indicator light on. If not press other keys within 15 seconds, automatically shut down, indicator light off.
5. Leave your skin 15 seconds, automatically shut down, indicator light off.
6. The power is on , press ON / INC button will begin to massage. The default vibration intensity is LEVEL 1, The default working mode is A, To select the other mode, press the PROGRAM button repeatedly according to the reciprocal diagram of indicator light and working mode.
7. When it reaches the set time or leave your skin to temporarily stop, must press ON/INC button to turn on the power. It returns to work, LED is bright but not blinking, Stick the massage patch on your skin, press ON/INC button repeatedly, it begin to work and indicator light blinks, working mode is the one before pausing, Continue pressing this button, the vibration mode will increase gradually and reach the highest vibration mode LEVEL 10 eventually.

OFF/DEC BUTTON MODE BUTTON

ON/INC BUTTON

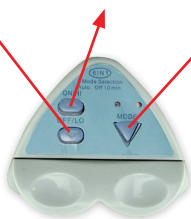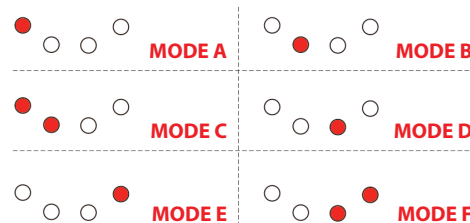RECIPROCAL DIAGRAM OF INDICATOR LIGHT AND WORKING MODE

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AB ENERGIZER GEBRAUCHSANWEISUNG